

## **Product Information Sheet**

# Anti-Phospho-CREB[Ser133]-176Yb

Catalog #: 3176005AClone: 87G3Package Size: 50 testsIsotype: Rabbit IgG


**Storage:** Store product at 4°C. Do not freeze. **Formulation:** Antibody stabilizer with 0.05% Sodium Azide

Cross Reactivity: Rat, Mouse, Human

### **Technical Information**

**Validation:** Each lot of conjugated antibody is quality control tested by CyTOF<sup>®</sup> analysis of stained cells using the appropriate postive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1  $\mu$ l for up to 3 X 10<sup>6</sup> live cells in 100  $\mu$ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

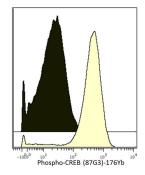

Human Jurkat T cells were incubated for 15 minutes in media alone (top) or with pervanadate (bottom). Cells were then fixed, permeabilized, and stained with 176Yb-anti-Phospho-CREB (87G3).

## Description

CREB is a 43 kDa leucine zipper transcription factor that belongs to the CREB/ATF family and regulates proliferation, differentiation and survival in diverse cell types including cells of neuronal and hematopoietic origin. The serine 133 residue within the kinase-inducible domain (KID) is phosphorylated by various kinases and this phosphorylation promotes the interaction of CREB with various transcriptional co-activators, in particular the histone acetyltransferases CREB-binding protein (CBP) or p300. Phosphorylation of CREB occurs in response to stimuli including growth factors, neurotransmitters and stress signals that increase intracellular cAMP or calcium levels. When activated, CREB binds to cAMP response elements, thereby initiating transcription of CREB target genes. The 87G3 antibody reacts with CREB phosphorylated at the serine 133 residue.

#### References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. Analytical Chemistry 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

#### **Contact Information:**

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

This product contains antibodies manufactured by and sold under license from Cell Signaling Technology and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

