

Anti-Human CD45 (HI30)-196Pt

Catalog number, package size: 3196001B, 100 tests
3196001C, 25 tests

Clone: HI30

Other Names: LCA, T200

Isotype: Mouse IgG1, κ

Reactivity: Human, Chimpanzee

Tag: 196Pt

Formulation: Antibody stabilizer with 0.05% sodium azide

Storage: Store at 2–8 °C. Do not freeze.


Application: CyTOF® suspension mass cytometry

Technical Information

Description: CD45, also known as leukocyte common antigen (LCA) and T200, is a type I transmembrane glycoprotein. It is expressed on the plasma membrane of all hematopoietic cells except mature red blood cells, platelets, and some plasma cells. Its intracellular domain is a tyrosine phosphatase that serves to regulate signal transduction in most hematopoietic cells.

Application: The platinum (Pt)-labeled antibody is designed and formulated for the application of suspension mass cytometry using the Fluidigm CyTOF suspension systems on healthy human PBMC.

Validation: Each lot of conjugated antibody is quality control-tested by CyTOF suspension mass cytometry analysis of stained cells using appropriate positive and negative cell staining and/or activation controls.

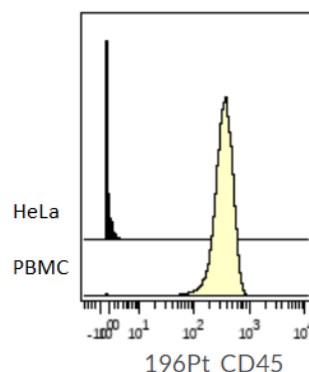

Human HeLa cells (top) and human PBMC (bottom) were stained with anti-CD45 (HI30)-196Pt. Total viable cells are displayed in the analysis.

Recommended use: Use 1 μ L for up to 3 million live cells in 100 μ L staining volume. We recommend titrating the antibody for optimal performance for each of the desired applications. Centrifuge the stock antibody at 12,000 \times g for 5 min to sediment antibody aggregates. Fixation is typically used in intracellular staining protocols or in barcoding with the Cell-ID™ 20-Plex Pd Barcoding Kit (201060). However, fixing before antibody staining can affect epitope structure and antibody binding, with the impact varying on the type and concentration of fixative and the protocol used. It is therefore important to perform a small, preliminary antibody staining experiment, with and without fixation, using noncritical samples.

NOTE This product is not compatible with Cell-ID Cisplatin and monoisotopic Cell-ID Cisplatin reagents of the same isotope. We strongly recommend that you titrate Cell-ID Intercalator-Ir for use with this product. You can use the Cadmium and Platinum Purity Matrices (FLDM-00038) or the Maxpar® Panel Designer tool to determine the metal purity or anticipated signal overlap, respectively.

Applicable Protocols

Before using this product, refer to the instructions in the Maxpar Cell Surface Staining with Fresh Fix Protocol (400276).

References

Bandura, D.R. et al. "Mass cytometry: technique for real time single cell multitarget immunoassay based on inductively coupled plasma time-of-flight mass spectrometry." *Analytical Chemistry* 81 (2009): 6,813–22.

Bendall, S.C. et al. "Single-cell mass cytometry of differential immune and drug responses across a human hematopoietic continuum." *Science* 332 (2011): 687–96.

Bodenmiller, B. et al. "Multiplexed mass cytometry profiling of cellular states perturbed by small-molecule regulators." *Nature Biotechnology* 30 (2012): 858–67.

Ornatsky, O.I. et al. "Highly multiparametric analysis by mass cytometry." *Journal of Immunological Methods* 361 (2010): 1–20.

Safety

Use standard laboratory safety protocols. Read and understand the safety data sheets (SDSs) before handling chemicals. To obtain SDSs, go to fluidigm.com and search for SDS-00025.

For technical support visit techsupport.fluidigm.com. | For general support visit fluidigm.com/support.

For Research Use Only. Not for use in diagnostic procedures.

Information in this publication is subject to change without notice. **Patent and license information:** fluidigm.com/legal/notices. **Limited Use Label License:** The purchase of this Fluidigm Instrument and/or Consumable product conveys to the purchaser the limited, non-transferable right to use with only Fluidigm Consumables and/or Instruments respectively except as approved in writing by Fluidigm. **Trademarks:** Fluidigm, the Fluidigm logo, Cell-ID, CyTOF, and Maxpar are trademarks and/or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2021 Fluidigm Corporation. All rights reserved. 05/2021