

Anti-Caspase3-172Yb

Catalog: 3172023A

Package Size: 50 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Rat, Mouse, Human, Monkey

Clone: 5A1E


Isotype: IgG

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3×10^6 live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

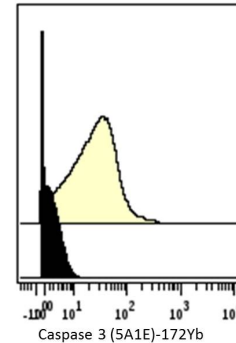

Human Jurkat T cells were incubated for 18 hours in media alone (bottom) or with etoposide (top). Cells were then fixed, permeabilized, and stained with 172Yb-anti-Cleaved Caspase 3 (5A1E).

Description

Caspase-3 (CPP32/YAMA/APOPAIN) is a 32 kDa cysteine protease that is activated during the early stages of apoptosis. Activation of caspase-3 requires proteolytic processing of its inactive zymogen into activated p17 and p12 fragments, which associate to form the active enzyme. Caspase-3 is either partially or totally responsible for the proteolytic cleavage of many key proteins, such as the nuclear enzyme poly (ADP-ribose) polymerase (PARP). Caspase-3 is widely distributed including high expression in cells of lymphoid origin, and active caspase-3 is a marker for cells undergoing apoptosis. Caspase-3 also plays an important role in morphogenetic cell death during development of the mammalian brain. The D3E9 antibody reacts with the active form of Caspase-3.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL **United States and other countries not in Europe or Asia** techsupport@fluidigm.com **Europe** techsupporteurope@fluidigm.com | **Asia** techsupportasia@fluidigm.com
 PHONE **United States (toll-free)** +1 866 358 4354 | **Europe** +33 1 60 92 42 40 | **Japan** +81 3 3662 2150 | **China** +86 21 3255 8368 | **All other countries** +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from Cell Signaling Technology and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015