

Anti-Human CD14-156Gd

Catalog: 3156019B Clone: HCD14
Package Size: 100 tests Isotype: IgG1


Storage: Store product at 4°C. Do not freeze. Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Human

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by $CyTOF^{\textcircled{R}}$ analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3 X 10⁶ live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

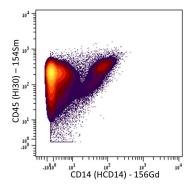

Human PBMCs stained with 154Sm-anti-CD45 (HI30) and 156Gd-anti-CD14 (HCD14).

Description

CD14, also known as the LPS co-receptor, is a glycosylphosphatidylinositol (GPI)-linked membrane glycoprotein. CD14 is expressed at high levels on antigen presenting innate immune cells including monocytes and macrophages, and at lower levels on granulocytes. It has also been reported to be expressed on some dendritic cell sub-populations

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly Multiparametric Analysis by Mass Cytometry. J Immunol Methods 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL United States and other countries not in Europe or Asia techsupport@fluidigm.com Europe techsupporteurope@fluidigm.com | Asia techsupportasia@fluidigm.com PHONE United States (toll-free) +1 866 358 4354 | Europe +33 1 60 92 42 40 | Japan +81 3 3662 2150 | China +86 21 3255 8368 | All other countries +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015