

Anti-Human CD16-165Ho

Catalog: 3165007B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Reactivity: Human

Clone: B73.1


Isotype: IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 μ l for up to 3×10^6 live cells in 100 μ l. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

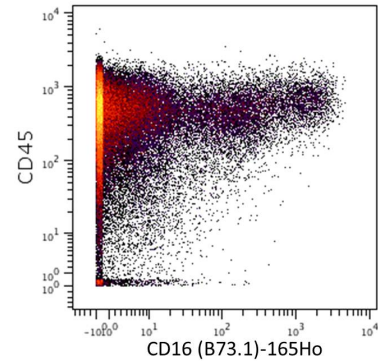

Human PBMCs stained with ¹⁵⁴Sm-anti-CD45 (HI30) and 165Ho-anti-CD16 (B73.1).

Description

CD16, also known as low affinity IgG receptor III (Fc γ RIII), is a polypeptide-anchored transmembrane protein. The CD16a isoform is expressed on the surface of NK cells, activated monocytes, macrophages, and placental trophoblasts. The CD16b isoform is expressed specifically on neutrophils. The B73.1 antibody binds to CD16-positive neutrophils with lower intensity when compared with some other CD16-specific antibodies.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

EMAIL **United States and other countries not in Europe or Asia** techsupport@fluidigm.com **Europe** techsupporteurope@fluidigm.com | **Asia** techsupportasia@fluidigm.com
 PHONE **United States (toll-free)** +1 866 358 4354 | **Europe** +33 1 60 92 42 40 | **Japan** +81 3 3662 2150 | **China** +86 21 3255 8368 | **All other countries** +1 650 266 6100

For Research Use Only. Not for use in diagnostic procedures.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 01/2015