

## Anti-Human CD194[CCR4]-175Lu

**Catalog #:** 3175021A

**Package Size:** 50 tests

**Storage:** Store product at 4°C. Do not freeze.

**Cross Reactivity:** Human

**Clone:** 205410


**Isotype:** Mouse IgG2b

**Formulation:** Antibody stabilizer with 0.05% Sodium Azide

### Technical Information

**Validation:** Each lot of conjugated antibody is quality control tested by CyTOF<sup>®</sup> analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1 µl for up to 3 X 10<sup>6</sup> live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

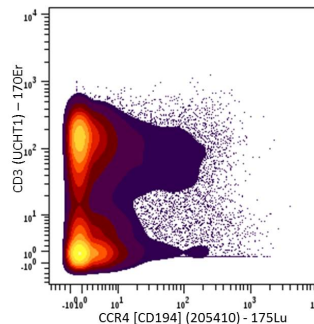

Human PBMCs stained with 170Er-anti-CD3 (UCHT1) and 175Lu-anti-CD194/CCR4 (205410).

### Description

CD194 also known as CCR4 (CC-chemokine receptor type 4) belongs to the G protein-coupled receptor family. It is expressed on stimulated and unstimulated T, and B cells, basophils, monocytes, and NK cells. Human peripheral blood Tregs can be divided into two distinct populations based on the expression of CCR4. Freshly isolated Tregs express CCR4 and presumably represent memory-type Tregs, and CCR4(-) Tregs require anti-CD3 Ab-mediated activation to acquire a regulatory activity. Depletion of CCR4(+) T cells leads to Th1-type polarization of CD4(+) T cells and augmentation of CD8(+) T cell responses to tumor Ags.

### References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

### Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com  
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

