

Anti-Human CD45RO-165Ho

Catalog #: 3165011B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Cynomolgus Monkey, Rhesus, Chimpanzee

Clone: UCHL1


Isotype: Mouse IgG2a

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3×10^6 live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

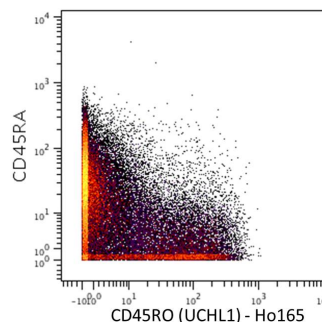

Human PBMCs stained with 153Eu-anti-CD45RA (HI100) and 165Ho-anti-CD45RO (UCHL1). Lymphocytes are displayed in the analysis.

Description

CD45RO, a splice variant of CD45, is a single chain membrane glycoprotein in the protein tyrosine phosphatase (PTP) family. It is expressed on activated memory T cells, B cell subsets, activated monocytes/macrophages and granulocytes.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Newell, E.W., et al. Cytometry by Time-of-Flight Shows Combinatorial Cytokine Expression and Virus-Specific Cell Niches within a Continuum of CD8+ T Cell Phenotypes. *Immunity* 36 January 2012: 142-152.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.
This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

