

Anti-Human CD71-175Lu

Catalog #: 3175011B

Package Size: 100 tests

Storage: Store product at 4°C. Do not freeze.

Cross Reactivity: Human

Clone: OKT-9


Isotype: Mouse IgG1

Formulation: Antibody stabilizer with 0.05% Sodium Azide

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

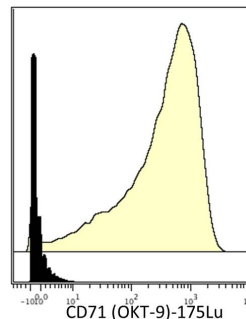

Human PBMCs were incubated for 3 days in media alone (bottom) or with PHA (top). Cells were then stained with 175Lu-anti-CD71 (OKT-9).

Description

CD71 is an integral membrane protein that mediates the uptake of transferrin-iron complexes. CD71 is a homodimeric glycoprotein that binds to diferric transferrin at the cell surface, and is most highly expressed on placental syncytiotrophoblasts, myocytes, basal keratinocytes, hepatocytes, endocrine pancreas, spermatocytes, and erythroid precursors. The maturation of erythroid precursors to erythrocytes results in loss of transferrin receptor expression, in concert with down-regulation of the machinery for hemoglobin synthesis. The high level of transferrin receptor within erythroid precursors makes it an excellent potential marker for evaluating erythroid elements within the bone marrow. CD71 has been studied for use in the assessment of leukemias and myelodysplastic syndromes (MDSs) in flow cytometric studies, in which expression serves as a lineage marker for the diagnosis of erythroid leukemia.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

Contact Information:

Sales: sales@DVSsciences.com | Support: support@DVSsciences.com
www.DVSsciences.com | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.
This product contains antibodies manufactured by and sold under license from eBioscience and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

