

Anti-Mouse CD117 (c-kit)-166Er

Catalog: 3166004B

Clone: 2B8

Package Size: 100 tests

Isotype: Rat IgG2b

Storage: Store product at 4°C. Do not freeze.


Formulation: Antibody stabilizer with 0.05% Sodium Azide

Reactivity: Mouse,

Technical Information

Validation: Each lot of conjugated antibody is quality control tested by CyTOF[®] analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

Recommended Usage: The suggested use is 1 µl for up to 3 X 10⁶ live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

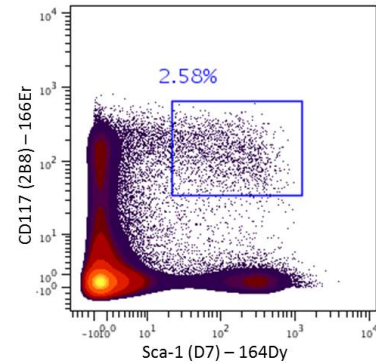

Mouse bone marrow cells stained with 164Dy anti-Sca-1 (D7) and 166Er anti-CD117 (2B8). Viable lineage- cells are displayed in the analysis. Lin-Sca-1+c-Kit+ LSK cells are indicated in the gate.

Description

CD117, also known as c-Kit and stem cell factor receptor, is an immunoglobulin superfamily member. It is a transmembrane receptor tyrosine kinase that binds stem cell factor (also known as steel factor and mast cell growth factor). CD117 is expressed on hematopoietic stem cells, progenitors committed to myeloid and/or erythroid lineages, T and B cell precursors, and mast cells.

References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. **Highly Multiparametric Analysis by Mass Cytometry.** *J Immunol Methods* 361 (1-2):1-20, 2010.

For technical support visit fluidigm.com/support

North America +1 650 266 6100 | Toll-free: +1 866 358 4354 in the US | support.northamerica@fluidigm.com **Europe** +33 1 60 92 42 40 | support.europe@fluidigm.com

China (excluding Hong Kong) +86 21 3255 8368 | techsupportchina@fluidigm.com **Japan** +81 3 3662 2150 | techsupportjapan@fluidigm.com

All other Asian countries +1 650 266 6100 | techsupportasia@fluidigm.com **Central and South America** +1 650 266 6100 | techsupportlatam@fluidigm.com

For Research Use Only. Not for use in diagnostic procedures.

This product contains antibodies manufactured by and sold under license from BioLegend and licensees thereof.

Information in this publication is subject to change without notice. **Safety data sheet information** fluidigm.com/sds **Patent and license information** fluidigm.com/legalnotices | Fluidigm, the Fluidigm logo, and CyTOF are trademarks or registered trademarks of Fluidigm Corporation in the United States and/or other countries. © 2015 Fluidigm Corporation. All rights reserved. 07/2015