

## Anti-Phospho-p44/42 MAPK (Erk1/2)[T202/Y204]-167Er

**Catalog #:** 3167005A

**Package Size:** 50 tests

**Storage:** Store product at 4°C. Do not freeze.

**Cross Reactivity:** Rat, Mouse, Human, Bovine, Canine, Porcine, Hamster, Monkey

**Clone:** D13.14.4E


**Isotype:** Rabbit IgG

**Formulation:** Antibody stabilizer with 0.05% Sodium Azide

### Technical Information

**Validation:** Each lot of conjugated antibody is quality control tested by CyTOF® analysis of stained cells using the appropriate positive and negative cell staining and/or activation controls.

**Recommended Usage:** The suggested use is 1 µl for up to 3 X 10<sup>6</sup> live cells in 100 µl. It is recommended that the antibody be titrated for optimal performance for each of the desired applications.

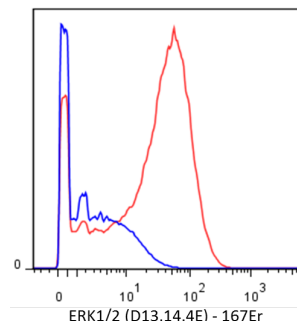

Human Jurkat T cells were incubated for 15 minutes with wortmannin (blue) or with PMA and ionomycin (red). Cells were then fixed, permeabilized, and stained with 167Er-anti-pERK1/2 (T202/Y204)(D13.14.4E).

### Description

ERK1 and ERK2, also known as p44 and p42 MAPKs, are similar (85% sequence identity) members of the mitogen activated protein kinase (MAPK) family of serine/threonine protein kinases. ERK1/2 signaling is important in the cellular response to a wide range of stimuli including growth factors, cytokines, and mitogens. The signal cascade upstream of ERK1/2 typically begins with receptor tyrosine kinases phosphorylating members of the Raf family and other MAP kinase kinase kinases (MAP3Ks), which thereby activate MEK1 and MEK2, the MAP kinase kinases (MAPKKs) directly responsible for phosphorylation of ERK1 and ERK2. ERK1 and ERK2 are activated through phosphorylation of the activation loop residues Thr202/Tyr204 and Thr185/Tyr187, and dual phosphorylation is required for full activity. ERK1/2 can activate the RSK family of kinases in the cytoplasm and several transcription factors (i.e. Elk-1) in the nucleus.

### References

Bandura, D. R., et al. Mass Cytometry: Technique for Real Time Single Cell Multitarget Immunoassay Based on Inductively Coupled Plasma Time-of-Flight Mass Spectrometry. *Analytical Chemistry* 81:6813-6822, 2009.

Ornatsky, O. I., et al. Highly multiparametric analysis by mass cytometry. *J Immunol Methods* 361 (1-2):1-20, 2010.

### Contact Information:

Sales: [sales@DVSsciences.com](mailto:sales@DVSsciences.com) | Support: [support@DVSsciences.com](mailto:support@DVSsciences.com)  
[www.DVSsciences.com](http://www.DVSsciences.com) | For assistance by phone: 855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.  
This product contains antibodies manufactured by and sold under license from Cell Signaling Technology and licensees thereof.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners. DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

