


Nebulizer Fit Kit

Catalog#: 101497
 Package Size: 1 unit
 Instrument: CyTOF and CyTOF 2

Technical Information

Description:

The Nebulizer Fit Kit works as pressure-reducing union connection for the nebulizer side arm and the nebulizer gas line through which Argon gas is supplied. The nylon Swagelok 1/4" nut connects to the gas line and the 3/16" nut holds the nebulizer side arm. The kit also contains a union body, a ferrule, a black O-ring and instructions for nebulizer connection.

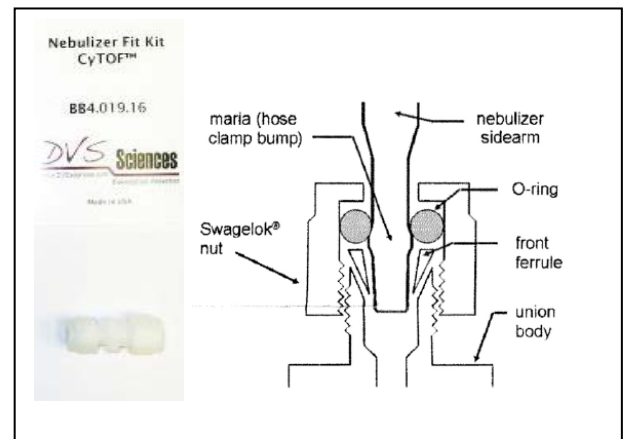

Contact Information:

Sales: sales@DVSsciences.com | Support: Support.NorthAmerica@DVSsciences.com or Support.Europe@DVSsciences.com
 www.DVSsciences.com | For assistance by phone: +1-855-DVS-CYTO

This product is intended for Research Use Only. This product is not intended for therapeutic or diagnostic purposes in humans or animals.

Vendor info & description quote as needed.

List supplier of raw material antibody here only if targets/clones are unique, licensed, and/or indicated as required in our agreement.

© 2011 DVS Sciences Inc. The trademarks mentioned herein are the property of DVS Sciences Inc. or their respective owners.

DVS Sciences Inc. products are warranted to meet stated product specifications and to conform to label descriptions when used and stored properly. Unless otherwise stated, this warranty is limited to 6 months from date of sales for products used, handled and stored according to DVS Sciences, Inc. instructions. DVS Sciences, Inc. sole liability for the product is limited to replacement of the product or refund of the purchase price. Product may not be resold, modified for resale or used to manufacture commercial products without prior written approval from DVS Sciences Inc.

